

# SOFTWARE DESIGN SPECIFICATIONS (SDS)

Project No : EMG-20-036

Date : 13.01.2021

Document Ref Number : EMG-20-036.3858.SDS

# AUTOMATIC ROTATIVE 6 HEAD (Pilfer Proof) CAPPING MACHINE

M-KAP-B01











EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 2

# **Document Approved**

## Prepared by

Name	Department/Function	Date	Sign

## Approved by

Name	Department/Function	Date	Sign

#### **Revised Date**

Rev.No.	Date	Prepared by	Description Changement
01	13.01.2021		Initial version







EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 3

1.	Purpose	4
2.	Scope	4
3.	System / Description of Equipment	4
3.1	Equipment / Description of System.	4
3.2	Machine Technical Datas.	4
4.	Description of Used System	5
4.1	System Architecture.	5
4.2	Software General View.	5
4.3	Programming Standard.	6
4.4	Software Informing	6
4.5	Programme Structure	7
4.6	Programming & Configuration and Software	7
4.6.1	TwinCAT 2	7
4.6.2	Movicon	8
5.	Access Authorization –User Management	9
6.	Changement Management	9
7.	Documentation.	9
8.	Comments and Explanation.	.10





#### Software Dizayn Spesifikasyonu (SDS)



EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 4

#### 1. Purpose

The purpose of this software design specification document is to define the design and function of the Rotative Plastering Capping Machine and its components project by the system designer

#### 2. Scope

This software design specification document is prepared for the Rotative Plastering Machine automation system and the automation system and its related components..

#### 3. System / Equipment Description

#### 3.1 Equipment / Description of System

The Rotative Capping Machine(Pilfer Proof) which is EMG-20-036 order code, consist of following equipment and components.

Equipment	Equipment	Equipment / Serial
Name	Code	No
Rotative Plastering Capping Machine	M-KAP-B01	3858

#### 3.2 Machine Technical Data

Machine	ROTATIVE CAPPING MACHINE
Туре	M-KAP-B01
Machine-No	3858
Weight	600 kg
Supply Voltage	3 Faz 380V 50Hz
Power Consumption	1.5 kW
Air Sply	4-6 bar








EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 5

#### 4. Description of Used System

### 4.1 System Architecture

Rotative Capping Machine

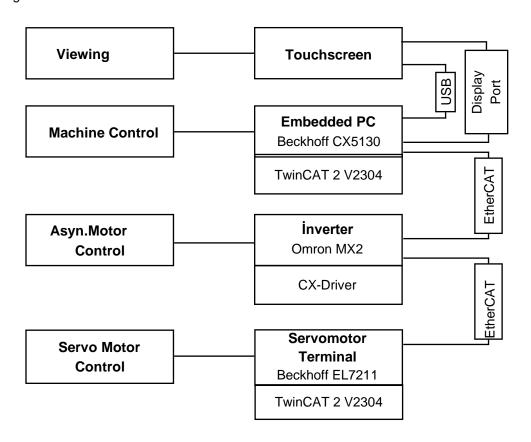

#### 4.2 Software General View

ROTATIVE CAPPING MACHINE		
Device	Embedded PC	
Manufacturer / Type	CX5130-0122	
Hardware Version	V3.8 2019-04-05	
Operating System	Windows Embedded Standart 32-Bit	
Operating System Version	6.1, Build 7601	





#### Software Dizayn Spesifikasyonu (SDS)



EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 6

#### 4.3 Programming Standard

The programming standard is based on the guidelines of GAMP 5 / ISO 9001.

#### 4.4 Software Informing

All control logic in the software is contained in the Beckhoff CX5130. Operator interface terminals read Embedded PC memory locations to monitor facility status and alarms. In addition, the operator terminal keeps all machine recipes / configuration parameters and saves them to the Embedded PC USB path or server path when necessary. All software modules and software sub-modules are called from the Embedded PC loop. The software modules used are a combination of program blocks and program function blocks.

#### **ROTATIVE CAPPING MACHINE**



Software type	TwinCAT / Movicon
Software name	M-KAP-B01
Software version	01
System request	RWin XP – Win 7 – Win 10
Data transfer	USB & Ethernet








EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 7

## 4.5 Programme Structure

### 4.6 Programming & Configuration and Software

#### 4.6.1 TwinCAT 2

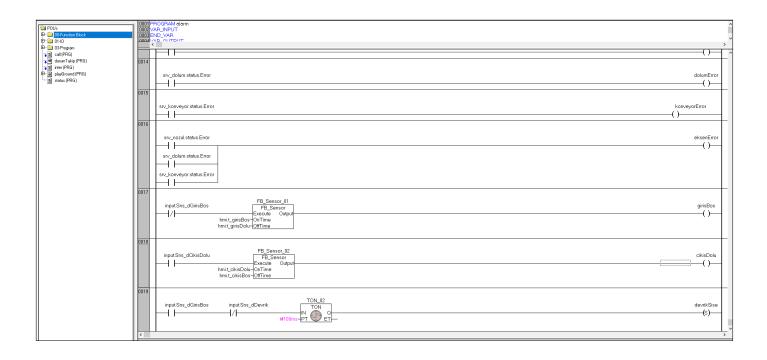

Description	
Software type	Programming Software
Software name	TwinCAT 2
Software version	2304
Job	TwinCAT programming software, build Embedded PC program and configure hardware is used.
System request	TwinCAT programming software works on a standard computer with Win XP, Win 7, Win10.
Data transfer	Application software data is transferred from the programming device to the Embedded PC via Ethercat.








EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 8

#### 4.6.2 Movicon

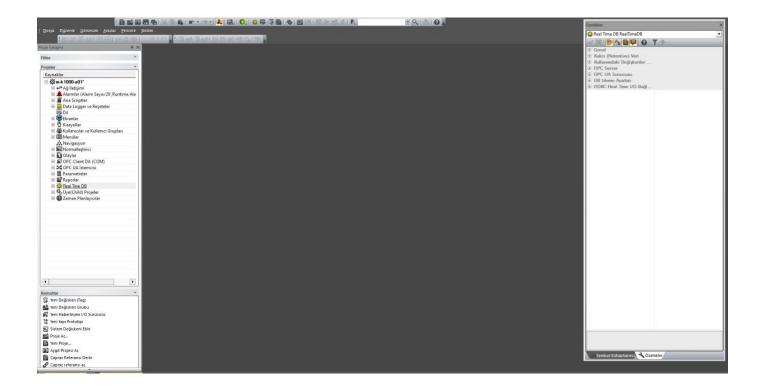

Description	
Software type	Display Software
Software name	Movicon 11.05
Software version	11.5.1185
Job	Movicon display software is used to visually prepare machine controls.
System request	Movicon screen software works on a standard computer with Win XP, Win 7, Win10.
Data transfer	Application software data is transferred from the programming device to the Embedded PC via Ethercat.





#### Software Dizayn Spesifikasyonu (SDS)



EMG-20-036.3858.SDS | Rev.No | 01 | Rev.Tarihi | 13.Ocak.2021 | Sayfa No | 9

#### 5. Access Authorization - User Management

Access authorization and user management are performed through groups and group Access Rights, which means that each user with their own password is assigned to a specific group and has access rights defined for that group. Some functions can be accessed without Group access rights, and the control system can be used unaffected. All other functions require logging in with a valid password. After logging in, all functions of the corresponding user group are available.

Automatic logoff occurs after 15 minutes.

User groups and passwords are managed through the SYSTEM /users menu. The following user groups with associated password are pre-set by ELEKTROMAG.

The following user groups and passwords are pre-set by ELEKTROMAG.

No	User group	Password
1	Keyless Input	NA
2	Operator	*****
3	Technician	*****
4	Administrator	****

Passwords will be notified by email to the company official.

#### 6. Changement Management

Change management will be handled in accordance with customer protocols and Quality Assurance System.

#### 7. Documentation

Program files for the machine'S PLC, HMI and frequency inverters are stored in the documentation folder marked with the data carrier symbol.

After the installation and validation of the entire system, it is recommended to the user to perform the "image import" process by the company.







EMG-20-036.3858.SDS Rev.No 01 Rev.Date 13.Jan.2021 Page No 10

# 8. Comments and Explanation

